


1 Definition of H.P

Helicobacter Pylori is a Gram-negative helicobacter that is parasitic on the human gastric mucosa.

Founded in 1983 by Warren and Marshall.

It can survive in a high degree of stomach acid, not only can coat stomach stick, but also pass through the protective barrier of stomach, causing stomach lesions.

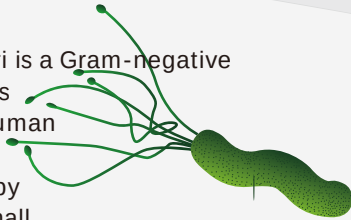

2 H.P infection do harm

The World Health Organization

Class I oncogenic source of gastric cancer



Detection rate of Gastric Ulcer is	Detection rate of Duodenal Ulcer is	Infection rate of Gastric Cancer is	The incidence of gastric cancer is increasing every year
70%	95%	89%	17.43%

3 Detection method of H.P

The breath test is the gold standard for detecting H.P

- Should be fasted for 4-6 hours before test. Take about 120ml warm water with 1 urea ¹⁴C capsule.
- Wait 10-20 minutes (Physical intened prohibited).
- Open the package of the breath card and the blowing nozzle, put the blow nozzle into the front section of the breath card.
- Blow into the breath card through blowing nozzle. The air should be blown out of the lung as smoothly as possible. you can change the air, but it is strictly prohibited to breathe in.
- Blow about 3-5 minutes, when the indicator window of the breath card changes from blue to white, stop blowing, if more than 3 minutes of incomplete color change, you can also stop blowing, at this time the breath card has absorbed saturation, does not affect the test results.
- Put the blowing nozzle into the trash can, give the breath card to the medical staff, and wait for the test results.

4 Route of infection

- mouth to mouth feeding
- fecal-oral transmission
- no separeate meal
- no serving chopsticks
- use same toothbrush
- kissing


5 H.P infection prevention

- 1 . Wash your hands before eating and after using the toilet.
- 2 . food is cooked at a high temperature.
- 3 . Eat less foods that irritate your stomach.
- 4 . Suggested separate meals.
- 5 . No mouth to mouth feeding.
- 6 . Replace toothbrush regularly.

轻轻松松吹口气检测 幽门螺旋杆菌

Blow it out easily
Detection for H. pylori

Our goal is to let the H.P nowhere to hide ! ! !



Henan Forever Medical Co.,Ltd



[¹⁴C] 呼气试验

¹⁴C Breath Test

Here is close to digestive health...